

OXYSURE™

With 24×7 Live
Oxygen Purity Display



Home 5 LPM
Oxygen Concentrator



Weight **13.5 kg**



Gold Standard Floating Ball Flowmeter

Provides precise oxygen flow adjustment in **0.1 LPM increments**, a trusted choice among healthcare professionals worldwide.

Suitable for Home & Clinical Use

Reliable Oxygen upto 5 Litres per minute

Delivers reliable **continuous oxygen flow from 0.5 to 5 LPM**, providing consistent therapy support across a wide range of prescribed requirements.



Built-in Nebulizer

Supports oxygen therapy and nebulization in one device, eliminating the need for a separate nebulizer.



One-Touch Therapy Start

Simply turn on the power switch and oxygen therapy begins instantly at your preset flow rate — no additional buttons required.




With Flow Memory During Power Cuts

Automatically restores the previously selected flow rate after power is restored, ensuring uninterrupted convenience.


Built for Continuous Operation

Designed for dependable, **long-duration oxygen therapy** to support everyday use whenever you need it.

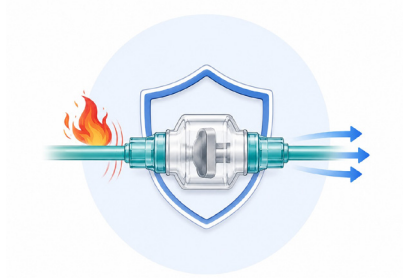

Fire Resisting Damper


Helps prevent the **spread of fire** through the oxygen delivery system, providing an additional layer of safety during oxygen therapy.

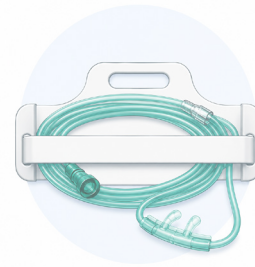

Cannula Holder Strap

Keeps the **nasal cannula** neatly secured when not in use, reducing tangles and improving organization.

Intuitive Control Panel

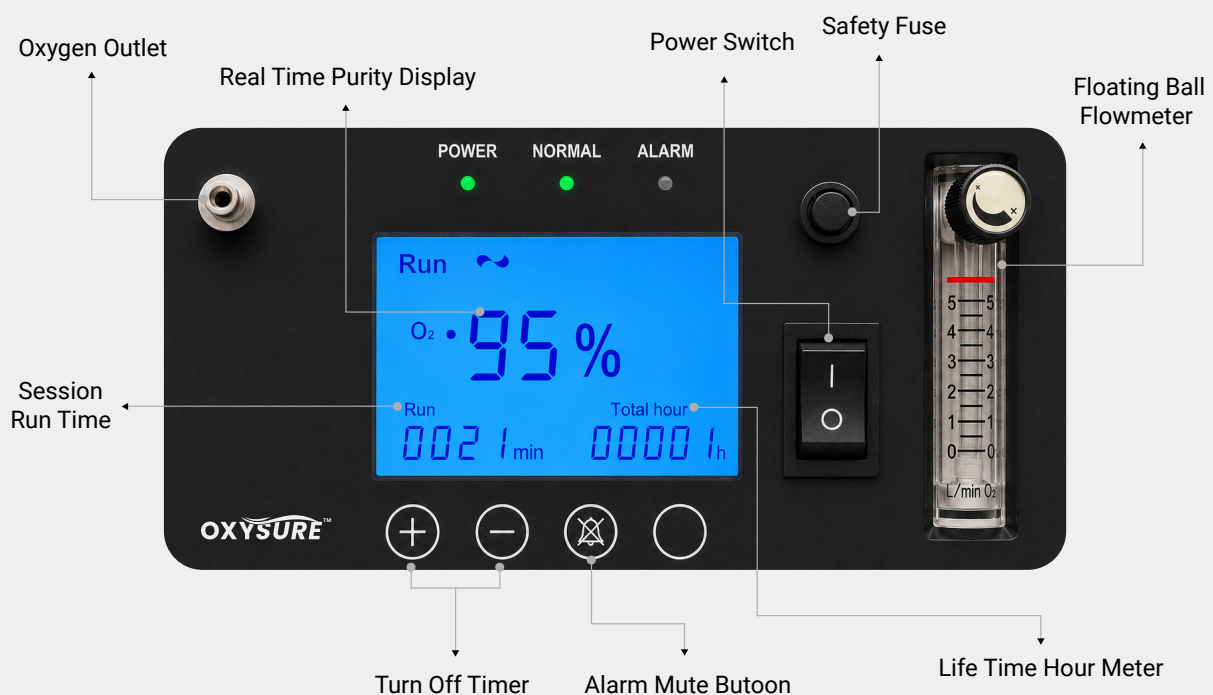

OXYSURE™ Home 5L

Reliable 5 LPM Oxygen, Wherever It's Needed.

Oxysure™ Home 5L is designed to deliver dependable continuous-flow oxygen for home care and clinical settings, with live oxygen purity monitoring and essential safety alarms.

Oxysure Home 5L Specifications

Model	Oxysure Home 5L
Oxygen Flow Rate	0.5–5 Litres Per Minute
Flowmeter	Floating Ball Type
Oxygen Concentration	93% ± 3%
Outlet Pressure	5.8–11.6 psi (40–80 kPa)
Power	320 W, 230 V AC, 50 Hz
Noise Level	40 dB (A)
Dimensions (L × W × H)	270 × 360 × 590 mm
Weight	13.5 kg
Conditions	Operating: 5°C to 40°C, ≤80% RH, Storage: -20°C to 55°C, ≤90% RH
Operating Altitude	Up to 7500 feet
Display Parameters	Session Run Time, Lifetime Hours, Purity Display, Turn-off Timer
Alarms	Power failure, System malfunction, and Low Oxygen Purity alarm
Nebulizer	Built-in port
Oxygen Purity Indicator	Green "Normal" indicator: Oxygen concentration ≥82% Red "Alarm" indicator: Oxygen concentration <82%
Oxygen Purity Display	Full Time Live Purity Percentage Display on screen



Oxysure 5L Oxygen Concentrator



Humidifier Bottle



Humidifier Connection Tube




7 feet Nasal Cannula



Nebulizer Mask & Tube



Firesafe Cannula Valve

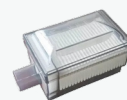

Filter Set Pre-installed



User Manual

Disclaimer: For oxygen therapy under the guidance of a qualified healthcare professional. For added safety, keep an adequate backup oxygen cylinder available for emergencies and unexpected interruptions in oxygen supply. Product specifications and features are subject to change without prior notice.

